

Multiplexed GuideRNA-expression to Efficiently Mutagenize Multiple Loci in *Arabidopsis* by CRISPR-Cas9

Julia Schumacher^{1,2}, Kerstin Kaufmann² and Wenhao Yan^{1,2,*}

¹Institute for Biochemistry and Biology, University of Potsdam, Potsdam, Germany; ²Department for Plant Cell and Molecular Biology, Institute for Biology, Humboldt-Universität zu Berlin, Berlin, Germany

*For correspondence: yanwen@uni-potsdam.de

[Abstract] Since the discovery of the CRISPR (clustered regularly interspaced short palindromic repeats)-associated protein (Cas) as an efficient tool for genome editing in plants (Li *et al.*, 2013; Shan *et al.*, 2013; Nekrasov *et al.*, 2013), a large variety of applications, such as gene knock-out, knock-in or transcriptional regulation, has been published. So far, the generation of multiple mutants in plants involved tedious crossing or mutagenesis followed by time-consuming screening of huge populations and the use of the Cas9-system appeared a promising method to overcome these issues. We designed a binary vector that combines both the coding sequence of the codon optimized *Streptococcus pyogenes* Cas9 nuclease under the control of the *Arabidopsis thaliana* *UBIQUITIN10* (*UBQ10*)-promoter and guideRNA (gRNA) expression cassettes driven by the *A. thaliana* *U6*-promoter for efficient multiplex editing in *Arabidopsis* (Yan *et al.*, 2016). Here, we describe a step-by-step protocol to cost-efficiently generate the binary vector containing multiple gRNAs and the Cas9 nuclease based on classic cloning procedure.

Keywords: CRISPR-Cas9, Multiplexing, gRNA, Gene editing, *Arabidopsis*

[Background] The RNA-guided Cas9-system is derived from the bacterial defense system against foreign DNA (Sorek *et al.*, 2013). It has been recognized as a method of choice for genome editing because of its high efficiency, easy handling and possibility of multiplex editing. In general, the Cas9-gene editing system involves a single synthetic RNA molecule, the gRNA that directs the Cas9 protein to target the desired DNA site for genome modification or transcriptional control. The gRNA-Cas9 complex recognizes the targeted DNA by gRNA-DNA pairing and requires the presence of a protospacer-adjacent motif (PAM). The PAM is represented by the nucleotides NGG or less specific NAG (with N for any nucleotide) in the target site following the gRNA-DNA pairing region. Thus, the approximate 20 nucleotide long gRNA spacer sequence, *i.e.*, the part of the gRNA sequence complementary to the DNA target site, determines the specificity of the complex. In this protocol, we describe the details to generate a binary vector that contains both the gRNA and the Cas9 coding sequence by classic cloning (Figure 1). As our vector system allows for subsequent addition of further gRNAs, it can be used for multiplex editing of the *Arabidopsis* genome and to obtain multiple, stably inherited alleles. Strong expression of the Cas9 protein and the gRNAs especially in proliferating tissues is achieved by the use of the *A. thaliana* *UBQ10*- and the *U6*-promoter, respectively. First mutations can be detected in the T1 generation, and T-DNA- and Cas9-free mutant plants may already be selected in

the T2 generation. Besides the selection of the gRNA and the construction of the plasmid, we give an overview of efficient genotyping methods required for detection of small or large deletions.

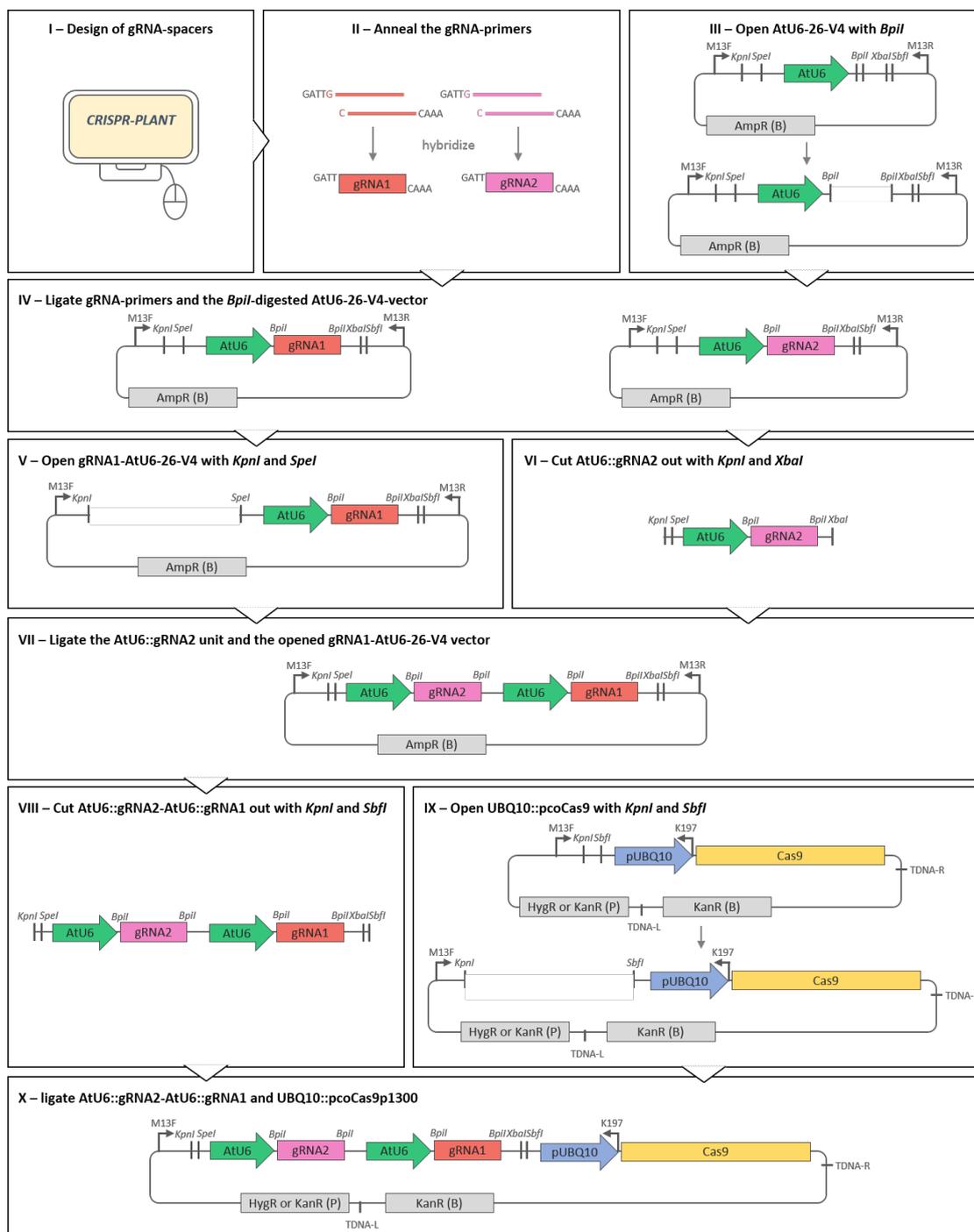


Figure 1. Scheme of the cloning procedure described in the protocol. Only a series of restriction and ligation steps is required to obtain a plant transformation vector equipped with a set of *AtU6*-driven gRNAs and the Cas9 enzyme under the control of the *UBQ10* promoter. Restriction enzyme cutting sites are displayed in italics. Arrows on vectors indicate primer

binding sites. TDNA-R and TDNA-L point out T-DNA right and left borders, respectively. (B) and (P) indicate selection markers for selection in bacteria or plants, respectively.

Materials and Reagents

Note: The protocol described here is based on the RNA-guided Cas9 system which was published recently (Yan et al., 2016). Only classical cloning methods such as restriction enzyme digestion and cohesive end ligation are required to construct plasmids ready for plant transformation.

1. Consumables

- a. 200 µl PCR tubes (Kisker Biotech, catalog number: G003-SF)
or 96-well PCR plates (SARSTEDT, catalog number: 72.1978.202)
StarSeal sealing tape (STARLAB, catalog number: E2796-9793)
- b. 1.5 ml reaction tubes (SARSTEDT, catalog number: 72.690.001)
- c. Pipette tips
- d. Scalpel
- e. Petri dishes, round, 9.2 x 1.6 cm (SARSTEDT, catalog number: 82.1472.001)
- f. Petri dishes, square, 10 x 10 x 2 cm (SARSTEDT, catalog number: 82.9923.422)
- g. Cuvettes for electroporation, e.g., Gene Pulser® cuvette 0.1 cm (Bio-Rad Laboratories, catalog number: 1652089)
- h. 5 ml glass pipette
- i. Micropore™ tape (VWR, catalog number: 115-8172)

Note: Any appropriate consumable can be used.

2. Competent cells

- a. One Shot® TOP10 chemically competent *Escherichia coli* (Thermo Fisher Scientific, Invitrogen™, catalog number: C4040)
- b. Electro-competent *Agrobacterium tumefaciens*, strain pGV3101 (for preparation of electro-competent *Agrobacteria*, please see Mersereau et al., 1990)

3. Plant material

- a. *Arabidopsis thaliana* Col-0

4. Plasmids

- a. AtU6-26-V4 (3.5 kb, ampicillin resistance marker [AmpR], available on request)
- b. UBQ10::pcoCas9p1300 (14.3 kb, kanamycin resistance marker [KanR] in bacteria, hygromycin or kanamycin resistance marker [HygR or KanR] in plants, available on request)
- c. Optional: pGEM®-T easy (Promega, catalog number: A3600)

5. Enzymes and buffers

Restriction

- a. *Bpi*I (Thermo Fisher Scientific, Thermo Scientific™, catalog number: ER1011)
- b. *Kpn*I-HF (New England Biolabs, catalog number: R3142)
- c. *Xba*I (New England Biolabs, catalog number: R0145)

- d. *SpeI*-HF (New England Biolabs, catalog number: R3133)
- e. *SbfI*-HF (New England Biolabs, catalog number: R3642)
- f. Cutsmart® buffer (New England Biolabs, catalog number: B7204S, supplied with the enzyme)

Ligation

- a. T4 DNA ligase (New England Biolabs, catalog number: M0202)
- b. T4 DNA ligase reaction buffer (New England Biolabs, catalog number: B0202S, supplied with the enzyme)

Genotyping (bacteria)

- a. Green *Taq* DNA polymerase (GenScript, catalog number: E00043)
 - b. 10x *Taq* buffer (GenScript, catalog number: B0005, supplied with the enzyme)
 - c. 10 mM dNTPs (Carl Roth, catalog number: K039)
 - d. Any DNA-loading dye, e.g., 6x gel loading dye, purple (New England Biolabs, catalog number: B7024)
 - e. Any DNA-size standard, e.g., GeneRuler™ 1 kb Plus DNA ladder (Thermo Fisher Scientific, Thermo Scientific™, catalog number: SM1331)
 - f. Agarose (Carl Roth, catalog number: 3810)
 - g. 10x TBE electrophoresis buffer (see Recipes)
 - Tris (Applichem, catalog number: A1379)
 - Boric acid (Applichem, catalog number: 131015)
 - EDTA (Applichem, catalog number: 131669)
6. Deionized water (sterile)
 7. Antibiotics (store stocks at -20 °C)
 - a. Ampicillin (Carl Roth, catalog number: K029, stock 50 mg/ml in deionized water)
 - b. Kanamycin (Carl Roth, catalog number: T832, stock 30 mg/ml in deionized water)
 - c. Rifampicin (Applichem, catalog number: A2220, stock 25 mg/ml in DMSO)
 - d. Gentamycin (Carl Roth, catalog number: 0233, stock 10 mg/ml in deionized water)
 - e. Hygromycin (Carl Roth, catalog number: CP13, stock 10 mg/ml in deionized water)
 8. Media (see Recipes)
 - a. YEB medium for *Agrobacterium*
 - Meat extract (Carl Roth, catalog number: 5770)
 - Yeast extract (Carl Roth, catalog number: 2904)
 - Peptone (Sigma-Aldrich, catalog number: 82303)
 - Sucrose (Carl Roth, catalog number: 4621)
 - Magnesium sulfate (MgSO₄) (Carl Roth, catalog number: P027)
 - Bacto-agar (Th. Geyer, CHEMSOLUTE®, catalog number: 9914-500G)
 - Sodium hydroxide (NaOH) (Merck, catalog number: 28245)
 - b. YT medium for *E. coli*
 - Sodium chloride (NaCl) (Carl Roth, catalog number: 9265)

- Yeast extract (Carl Roth, catalog number: 2904)
- Peptone (Sigma-Aldrich, catalog number: 82303)
- Bacto-agar (Th. Geyer, CHEMSOLUTE®, catalog number: 9914-500G)
- Sodium hydroxide (NaOH) (Merck, catalog number: 28245)
- c. ½ MS for *Arabidopsis*
 - MS + B5 Vitamins (Duchefa Biochemie, catalog number: M0231)
 - Potassium hydroxide (KOH) (Merck Millipore, catalog number: 105033)
- 9. DNA purification
 - a. Gel-purification, e.g., NucleoSpin® Gel and PCR clean-up (MACHEREY-NAGEL, catalog number: 740609)
 - b. Plasmid isolation, e.g., NucleoSpin® Plasmid EasyPure (MACHEREY-NAGEL, catalog number: 740727)
- 10. Silwet L-77 (Lehle seeds, catalog number: VIS-30) for plant transformation
- 11. Bleach (Carl Roth, catalog number: 9062)
- 12. 37% HCl (Merck Millipore, catalog number: 100317)
- 13. Genotyping (plants)
 - a. Optional: Phire Plant Direct PCR Kit (Thermo Fisher Scientific, Thermo Scientific™, catalog number: F-130WH)
 - b. Optional: T7 endonuclease I (New England Biolabs, catalog number: M0302) plus NEBuffer 2 (New England Biolabs, catalog number: B7002), 250 mM EDTA
 - c. Optional: pGEM®-T easy (Promega, catalog number: A3600)
- 14. Oligonucleotides (5'-3'), 10 pmol/μl
 - a. M13F_TGTAACGACGGCCAGT
 - b. M13R_CAGGAAACAGCTATGACC
 - c. K197_CTGTTAATCAGAAAACTCAG

Equipment

Note: No specific equipment is required. Any appropriate device can be used.

1. Computer with internet access
2. Thermocycler (e.g., Eppendorf, model: Mastercycler® nexus)
3. Thermoblock (e.g., Eppendorf, model: Thermomixer® comfort)
4. Plate incubators (28 °C, 37 °C)
5. Shakers (28 °C, 37 °C)
6. Horizontal gel-electrophoresis system (e.g., Bio-Rad Laboratories, model: Mini-Sub® Cell GT System)
7. UV transilluminator (e.g., Bio-Rad Laboratories, model: GelDoc™ XR+ System)
8. Electroporator, (e.g., Bio-Rad Laboratories, model: MicroPulser™ Electroporator)
9. 10 L desiccator

10. Fume hood
11. 200 ml beaker

Procedure

1. Preparation of gRNA-spacers with CRISPR-PLANT[®] (Xie *et al.*, 2014) and primer design (Figure 1I)

- a. Go to <http://www.genome.arizona.edu/crispr/index.html>.
- b. Select 'Search'.
- c. Select 'Select Species' and choose '*Arabidopsis thaliana*'.
- d. Select 'Chromosome' and enter the coordinates of the target region in 'From' and 'To'.
- e. Confirm with 'Search by region'.

You will be returned Class 0.0 and Class 1.0 gRNAs. Class 0.0 gRNAs are supposed to be more specific. The SeqID is the coordinate of the gRNA in your predefined region. The minMM_GG (and minMM_AG) value refers to the minimum amount of mismatches the 12 bp sequence adjacent to the PAM motif has to any other sequence in the *A. thaliana* genome. The higher the value the smaller is the risk for off-targeting. Therefore, sequences with a minMM_GG higher than three should be considered.

- f. Copy the optimal 20 bp long spacer sequence. Here, X represents any nucleotide.

Forward sequence: 5'-XXXXXXXXXXXXXXXXXXXX-3'

Reverse sequence: 5'-XXXXXXXXXXXXXXXXXXXX-3'

As the *U6*-promoter drives the expression of the gRNA, the initiation of transcription requires a guanine (G) as the transcription start site. If the chosen spacer sequence does not start with a 'G', change the first base. This will have no considerable effects on the specificity of Cas9 as the specificity largely depends on the 12 bp adjacent to the PAM motif.

Forward sequence: 5'-GXXXXXXXXXXXXXXXXXXXX-3'

Reverse sequence: 5'-XXXXXXXXXXXXXXXXXXXXC-3'

- g. Add 'GATT' and 'AAAC' to the 5'-ends of the forward and reverse sequences, respectively, to generate *Bpi*I-compatible overhangs.

Forward primer: 5'-GATTGXXXXXXXXXXXXXXXXXXXX-3'

Reverse primer: 5'-AAACXXXXXXXXXXXXXXXXXXXXC-3'

- h. Now you have a complementary primer pair ready to order.

2. Anneal the gRNA-primers (Figure 1II)

- a. Mix 10 μ l of each primer (10 pmol/ μ l). Denature primers for 3 min at 98 °C. Hybridize primers in a thermocycler by decreasing the temperature from 98 °C to 22 °C. Pause at each degree for 30 sec.

3. Digest the AtU6-26-V4-vector with *Bpi*I according to manufacturer's instructions (Figure 1III)

- a. Mix the digest with DNA-loading dye and load it on a 1% agarose gel. Separate digested from undigested bands by gel electrophoresis.

- b. Purify the linear vector backbone (3,517 bp) from gel.
4. Ligate the annealed gRNA-primers and the *Bpi*I-digested AtU6-26-V4-vector in a 10 μ l ligation reaction (Figure 1IV).
 - a. Mix 25 ng to 50 ng vector, 1 μ l of 10x T4 DNA ligase reaction buffer and 0.5 μ l T4 DNA ligase (200 U). Fill the rest of the volume with the hybridized primer.
 - b. Incubate the reaction for 10 min at room temperature.

Note: Use 4 °C over night incubation if you encounter problems getting positive transformants in step 6.
5. Use 3 μ l of the ligation reaction to transform chemically competent TOP10 *E. coli* cells.
 - a. Thaw 50 μ l competent cells on ice.
 - b. Add 3 μ l of the ligation reaction to the cells.
 - c. Incubate on ice for 15 min.
 - d. Heat shock the cells for 45 sec at 42 °C in a thermoblock and place back on ice for 2 min.
 - e. Add 400 μ l of YT medium and incubate the cells in a thermoblock for 60 min at 37 °C and 750 rpm.
 - f. Plate 60 to 80 μ l of the transformation on solid YT medium with 50 μ g/ml ampicillin (YT-amp).
6. Do a colony-PCR to genotype transformants. Usually testing 5 to 10 colonies is sufficient to obtain positive transformants.
 - a. Set up a 20 μ l PCR reaction per colony as follows. Final concentrations are indicated in brackets:
 - 2 μ l 10x *Taq* buffer (1x)
 - 0.6 μ l M13F (0.5 μ M)
 - 0.6 μ l gRNA-specific reverse primer (0.5 μ M)
 - 0.4 μ l dNTPs (0.2 mM)
 - 0.1 μ l Green *Taq* DNA polymerase (0.5 U)
 - 16.3 μ l deionized water
 - b. Use a pipette tip to pick a colony and restreak the cells on an YT-amp plate. Directly afterwards dip the tip into the PCR reaction to transfer a small amount of cells. Incubate the plate at 37 °C for approximately 5 h and run the PCR in the thermocycler using the following conditions:

	94 °C	hold	
Initial Denaturation	94 °C	3.15 min	
Denaturation	94 °C	30 s	} 30 cycles
Annealing	54-60 °C	45 s	
Elongation	72 °C	35 s	
Final Extension	72 °C	7 min	
	4 °C	hold	

Note: The annealing temperature depends on the gRNA-specific reverse primer. M13F works well within the range of 54 °C to 60 °C during annealing. Instead of M13F and the

gRNA-specific reverse primer you can also use M13R and the gRNA-specific forward primer with an elongation time of 20 sec.

- c. Add 4 μ l of 6x DNA loading dye to the reaction and load it on a 1% agarose gel. The expected amplicon size is 581 bp for the PCR reaction using M13F and the gRNA-specific reverse primer or 257 bp for the reaction with M13R and the gRNA-specific forward primer.
 - d. Prepare 3 ml of YT-amp and pick the PCR-positive clones from the restreak plate. Incubate over night at 37 °C with 190 rpm.
 - e. Isolate the plasmid from positive clones and check for the correct insert by Sanger sequencing using M13R as sequencing primer.
7. Optional (for multiplexing):
- If more than one gRNA is needed, prepare each gRNA in a separate AtU6-26-V4-vector first (gRNA1-AtU6-26-V4, gRNA2-AtU6-26-V4, etc.) and combine them into one single gRNA2-gRNA1-AtU6-26-V4-vector.
- a. Cut gRNA1-AtU6-26-V4 open with *KpnI*-HF and *SpeI*-HF (3,519 bp, Figure 1V).
 - b. Cut gRNA2 out of gRNA2-AtU6-26-V4 using *KpnI*-HF and *XbaI* (643 bp, Figure 1VI).
 - c. Load reactions on gel, gel purify positive bands.
 - d. Ligate *KpnI*-gRNA2-*XbaI* and *SpeI*-gRNA1-AtU6-26-V4-*KpnI* (*XbaI* and *SpeI* produce compatible ends, Figure 1VII).
 - e. Transform TOP10 *E. coli* cells, select transformants on YT-amp plates.
 - f. Genotype clones with M13F and a specific reverse primer for the last added gRNA (amplicon size = 581 bp).
 - g. Isolate plasmids and sequence with M13R.
Repeat this step to add more gRNAs.
8. Digest gRNA2-gRNA1-AtU6-26-V4 with *KpnI*-HF and *SbfI*-HF according to manufacturer's instructions. With this step, you will cut out the gRNA2-gRNA1 unit including U6-promoters (Figure 1VIII).
- a. Mix the digest with DNA-loading dye and load it on a 1% agarose gel. Separate digested from undigested bands by gel electrophoresis.
 - b. Purify the gRNA-fragment (660 bp for a single gRNA-fragment, 1,287 bp for a double gRNA-fragment) from gel.
9. Digest UBQ10::pcoCas9p1300 with *KpnI*-HF and *SbfI*-HF according to manufacturer's instructions to open the vector (Figure 1IX).
- a. Mix the digest with DNA-loading dye and load it onto a 1% agarose gel. Separate digested from undigested bands by gel electrophoresis.
 - b. Purify the linearized vector backbone (14.3 kb) from gel.
10. Ligate the gRNA2-gRNA1 unit and UBQ10::pcoCas9p1300 in a 20 μ l reaction (Figure 1X).
- a. Mix 25 ng to 50 ng vector, 10 ng to 20 ng gRNA2-gRNA1, respectively, 2 μ l of 10x T4 DNA ligase reaction buffer, deionized water ad 19 μ l, and 1 μ l T4 DNA ligase (400 U).
 - b. Incubate the reaction at 4 °C over night.

11. Use 3 μ l of the ligation reaction to transform chemically competent TOP10 *E. coli* cells. Select transformants on YT plates with 30 μ g/ml kanamycin.
12. Do a colony-PCR to genotype transformants as described in step 6.
 - a. Genotype clones with K197 and a specific forward primer for the first added gRNA (here gRNA1, product size 801 bp).
 - b. Isolate and sequence plasmids with M13F (optional with K197).
13. Use 50 ng AtU6::gRNA2-AtU6::gRNA1-UBQ10::pcoCas9p1300 to transform *A. tumefaciens* via electroporation.
 - a. Thaw 50 μ l competent *A. tumefaciens* on ice.
 - b. Pre-cool the cuvette on ice.
 - c. Add 50 ng final vector to the *A. tumefaciens* cells, transfer to the cuvette and place the cuvette in the electroporator.
 - d. Choose the 'Agr' settings at the MicroPulser™ electroporator (Bio-Rad, 1 pulse at 2.2 kV) or equivalent settings to transform the cells.
 - e. Add 1 ml of YEB medium and transfer into a fresh 1.5 ml reaction tube. Incubate at 28 °C and 750 rpm for 2 to 3 h.
 - f. Plate 60 μ l of the transformation and select on YEB medium with 30 μ g/ml kanamycin, 100 μ g/ml rifampicin and 25 μ g/ml gentamycin at 28 °C for 2 days.
14. Use the transformed *A. tumefaciens* to transform *A. thaliana* by floral dip with infiltration (Clough and Bent, 1998; Chen, 2011)
15. Sterilize seeds, e.g., by vapor-phase sterilization.
 - a. Transfer seeds to 1.5 ml reaction tubes.
 - b. Place open tubes in a 10 L desiccator under the fume hood.
 - c. Place a 200 ml beaker filled with 100 ml bleach (6%) into the desiccator.
 - d. Use a glass pipette to add 3 ml HCl (37%) to the bleach and close desiccator lid quickly.
Note: Careful: The vapor can chemically burn skin and eyes!
 - e. Sterilize for 3 to 4 h.
16. Sow seeds on ½ MS with 20 μ g/ml hygromycin or kanamycin (depending on the resistance marker from the UBQ10::pcoCas9p1300 of your choice) to select for positive transformants. Seal plates with Micropore™ tape. Stratify seeds in darkness at 4 °C for 3 days. Germinate seeds in long day conditions (22 °C, 16:8 h light:dark photoperiod) until transformants reach 4 true-leaf stage. Transplant seedlings to soil and continue growth in long day conditions.
17. Genotype the transformants.
As in the T1 generation the activity and efficiency of the Cas9 enzyme but also of the DNA repair mechanisms might differ from cell to cell, each plant most likely will be a mosaic of different genotypes. Therefore we suggest to pool genotyping samples of one plant, such as pieces of cauline leaves from different branches. If a plant with the expected mutation was identified in T1, label and genotype the single branches as they can have different genotypes (Figures 2C and 2D). The method of genotyping depends on the anticipated results.

- a. Big fragment deletions which had been generated by multiple gRNAs, can be identified by PCR followed by gel electrophoresis of the PCR products (Figures 2A and 2B).

Note: Design primers flanking the expected deletion in a distance of approximately 200 bp to 400 bp to each gRNA target site. For high throughput genotyping we suggest to prepare 10 μ l PCR reactions using the Phire Plant Direct PCR Kit.

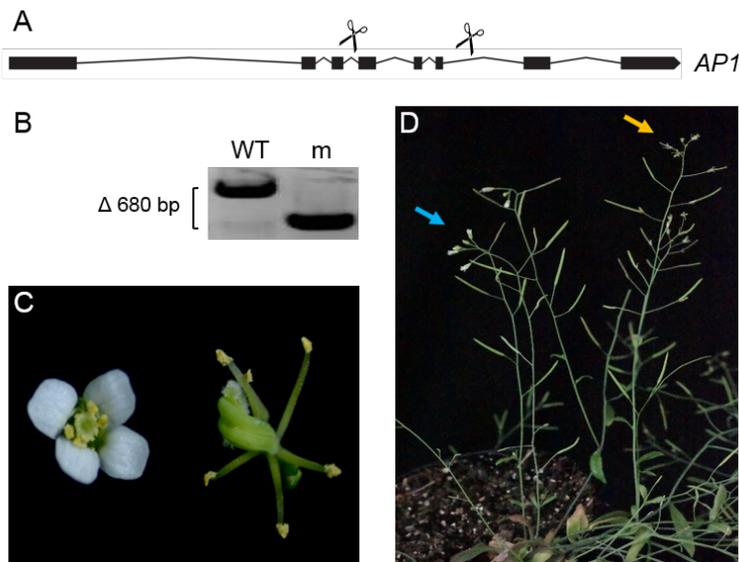


Figure 2. Detection of big fragment deletions in the *AP1* locus. A. The *AP1* locus was targeted using two different gRNAs (scissors). Target sites are sequences in the third and sixth intron of *AP1*. B. Genotyping of a phenotypically wild type (WT) and an *ap1*-like (m) branch in the T1 generation revealed the presence of a 680 bp deletion in the *AP1* locus in the *ap1*-like branch (m). C and D. Plants in the T1 generation can be a mosaic of different genotypes. Two different flower phenotypes, wild type (C, left; D, blue arrow) and *ap1* (C, right; D, yellow arrow) can be found on two different branches of the same plant.

- b. Small deletions or polymorphisms that produce a minimum mismatch of 2 bp to the wild type allele can be detected by T7 endonuclease digestion (Vouillot *et al.*, 2015).
- i. Perform the genotyping PCR as mentioned for big fragment deletions.
 - ii. Set up the hybridization reaction by adding 1.12 μ l NEBuffer 2 to the PCR product without purification of the DNA. Prepare WT control and sample separately. When pooling different branches in T1, it is not necessary to add wild type PCR product to the sample as we usually do not expect Cas9 efficiencies to be higher than the T7 endonuclease detection limit (Vouillot *et al.*, 2015). But if you already see a clear phenotype when genotyping single branches, you can replace 20% of the sample volume by wild type to not miss homozygous mutant branches.
 - iii. Hybridize PCR products using the following conditions:

95 °C	hold	
95 °C	5 min	
95 °C	1 s	} -2 °C/s
85 °C	1 s	
25 °C	1 s	} -0.1 °C/s
4 °C	hold	

- iv. Add 0.1 µl of T7 endonuclease (1 U) and incubate for 30 min at 37 °C.
 - v. Stop the reaction with 0.8 µl 250 mM EDTA.
 - vi. Separate cleavage products on a 2% agarose gel.
- c. Point mutations are not detected by T7 endonuclease digest efficiently (Vouillot *et al.*, 2015) and the IDT® Surveyor® nuclease (Integrated DNA Technologies), which is supposed to cut at point mutations, appeared to be incompatible with the Phire Plant Direct PCR Mix. Considering the economic factors of high throughput mutant screening, we did not increase the amount of Surveyor® enzyme beyond manufacturer’s recommendations. Here, we suggest to use the T7 endonuclease assay to exclude plants with mutations of 2 bp or more. Then PCR products from the rest of the samples can be purified and sent for Sanger sequencing.
18. Confirm the mutation by sequencing.
 19. Harvest the seeds from the branch that has a mutation.
 20. Sow the seeds to generate a small T2 population of approximately 100 individuals. Select the plant that carries the mutation but lost the T-DNA.
 - a. To identify plants that still carry the transgene, use M13F and the gRNA-specific reverse primer (product size 557 bp). Exclude PCR-positive plants.
 - b. Genotype for the expected mutation as described in step 17.

Note: To verify small fragment deletions in T2 (Figure 3), prepare WT, sample and a 1:1 mix of WT and sample prior to hybridization. The WT-sample-mix is necessary as plants can be homozygous for the mutant allele. In that case they do not produce heteroduplex DNA when hybridized and will not be digested by T7 endonuclease.

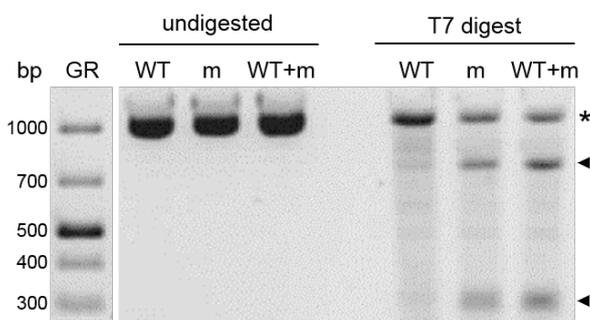


Figure 3. Detection of small deletions or polymorphisms by T7 digest. PCR samples of T2 wild type (WT), mutant (m) and a 1:1 mix of wild type and mutant (WT+m) were hybridized. Half

of the reaction volume was incubated without T7 endonuclease (undigested) and the other half was incubated with 1 U of the enzyme (T7 digest) and separated by gel electrophoresis on 2% agarose. The presence of cleavage products in the mutant sample is the result of at least two different alleles, which can be wild type and mutant or different mutant alleles or a mix of both options. As the addition of wild type to the mutant sample (WT+m) does not dilute the signal of the cleavage products, there are most likely different mutant alleles, which later was confirmed by sequencing (not shown). Asterisk indicates undigested PCR products, black arrows point out cleavage products. GR indicates the GeneRuler™ 1 kb Plus DNA ladder.

Note: Once the protocol is set up, the 'undigested' controls can be omitted.

Data analysis

Mutations obtained in the T1 generation can be confirmed by Sanger sequencing as mentioned in step 17. But, as neighboring cells might harbor different mutations and as the mutation might be biallelic, we recommend to gel-purify the mutant PCR product and subclone it to pGEM®-T or any other cloning vectors. Sequence five to ten positive colonies to get an overview of gained mutations. After the selection of a T-DNA-free mutant plant in the T2 generation (step 19) confirm the presence and type of the mutation again by sequencing.

Note: PCR products generated by the Phire Plant Direct PCR Kit are blunt ended and require A-tailing before ligation with pGEM®-T.

Notes

This protocol is simple and straightforward. So far, more than 50 binary vectors containing different sets of gRNAs have been generated by our group and used for genome editing of *Arabidopsis*.

Recipes

1. 10x TBE electrophoresis buffer (1 L)
 - 1 M Tris
 - 1 M boric acid
 - 0.02 M EDTA
 - Deionized water ad 1,000 ml
 - Dilute 10x TBE buffer with deionized water to 1:10 to use as gel running buffer
2. YEB medium for *Agrobacterium* (1 L)
 - 5 g beef extract
 - 1 g yeast extract
 - 5 g peptone
 - 5 g sucrose

- 500 mg MgSO₄
10 g Bacto-agar
Deionized water ad 1,000 ml
Adjust pH to 7.0 using NaOH
Autoclave YEB medium at 121 °C and 2 bar for 15 min
3. YT medium for *E. coli* (1 L)
5 g NaCl
5 g yeast extract
8 g peptone
(15 g Bacto-agar)
Deionized water ad 1,000 ml
Adjust pH to 7.0 using NaOH
Autoclave YT medium at 121 °C and 2 bar for 15 min
4. ½ MS for *Arabidopsis* (1 L)
2.3 g MS + B5 Vitamins
8 g Phyto-agar
Deionized water ad 1,000 ml
Adjust pH to 5.7 using KOH
Autoclave ½ MS at 121 °C and 2 bar for 15 min
Prepare sterile 50% sucrose (filter sterilize or autoclave separately)
Add 1% sucrose to ½ MS prior to pouring the plates
Note: The addition of sucrose is not required. The advantage is that seedlings recover better on selective media. The disadvantage is that it promotes fungal growth on the medium.

Acknowledgments

This protocol is based on our work previously published in Plant Methods ([Yan et al., 2016](#)). Work of J.S. was supported by the IMPRS-PMPG fellowship. K.K. wishes to thank the Alexander-von-Humboldt foundation and the BMBF for support. The authors declare that they have no competing interests.

References

1. Chen, L. (2011). [Simplified Agrobacterium-mediated transformation of Arabidopsis thaliana](#). *Bio-protocol* Bio101: e41.
2. Clough, S. J. and Bent, A. F. (1998). [Floral dip: a simplified method for Agrobacterium-mediated transformation of Arabidopsis thaliana](#). *Plant J* 16(6): 735-743.

3. Li, J. F., Norville, J. E., Aach, J., McCormack, M., Zhang, D., Bush, J., Church, G. M. and Sheen, J. (2013). [Multiplex and homologous recombination-mediated genome editing in *Arabidopsis* and *Nicotiana benthamiana* using guide RNA and Cas9](#). *Nat Biotechnol* 31(8): 688-691.
4. Mersereau, M., Pazour, G. J. and Das, A. (1990). [Efficient transformation of *Agrobacterium tumefaciens* by electroporation](#). *Gene* 90(1):149-51.
5. Nekrasov, V., Staskawicz, B., Weigel, D., Jones, J. D. and Kamoun, S. (2013). [Targeted mutagenesis in the model plant *Nicotiana benthamiana* using Cas9 RNA-guided endonuclease](#). *Nat Biotechnol* 31(8): 691-693.
6. Shan, Q., Wang, Y., Li, J., Zhang, Y., Chen, K., Liang, Z., Zhang, K., Liu, J., Xi, J. J., Qiu, J. L. and Gao, C. (2013). [Targeted genome modification of crop plants using a CRISPR-Cas system](#). *Nat Biotechnol* 31(8): 686-688.
7. Sorek, R., Lawrence, C. M. and Wiedenheft, B. (2013). [CRISPR-mediated adaptive immune systems in bacteria and archaea](#). *Annu Rev Biochem* 82: 237-266.
8. Vouillot, L., Thelie, A. and Pollet, N. (2015). [Comparison of T7E1 and surveyor mismatch cleavage assays to detect mutations triggered by engineered nucleases](#). *G3 (Bethesda)* 5(3): 407-415.
9. Xie, K., Zhang, J. and Yang, Y. (2014). [Genome-wide prediction of highly specific guide RNA spacers for CRISPR-Cas9-mediated genome editing in model plants and major crops](#). *Mol Plant* 7(5): 923-926.
10. Yan, W., Chen, D. and Kaufmann, K. (2016). [Efficient multiplex mutagenesis by RNA-guided Cas9 and its use in the characterization of regulatory elements in the *AGAMOUS* gene](#). *Plant Methods* 12: 23.