

Gliding Assay to Analyze Microtubule-based Motor Protein Dynamics

Albert Shim, Gin Tezuka, Luke Kupcha and Li Tao*

Department of Biology, University of Hawaii at Hilo, 200 West Kawili Street, Hilo, Hawaii 96720, USA

*For correspondence: ltao2@hawaii.edu

[Abstract] The purpose of this protocol is to provide an updated method of performing microtubule gliding assays and visualizing it using fluorescence microscopy.

Keywords: Gliding assay, Mitotic motor, *in vitro* motility, Polarity-marked microtubules, Glass chamber

[Background] Mitotic spindles are protein machinery that dominate mitosis. The mitotic spindle utilizes microtubule-based motor proteins to organize itself, and exert forces to drive cell division. Microtubule-based motor proteins produce mechanical work using energy derived from ATP hydrolysis (Coppin *et al.*, 1997). Motor proteins translocate microtubules in a unidirectional manner. The behavior of motility can be observed by *in vitro* gliding assay (Tao and Scholey, 2010), in which the motors are affixed onto a glass surface and supplied with microtubules and ATP. The motility of the motor proteins can then be studied using fluorescence microscopy and the details of their dynamic behavior can be observed in real time. This updated protocol will allow analysis of microtubule-based motor protein function with the use of *in vitro* microtubule gliding assays (Tao *et al.*, 2006 and 2016).

Materials and Reagents

1. Pipette tips (Corning, catalog number: 4860)
2. Centrifuge tubes (Corning, catalog number: 430290)
3. Thickwall polycarbonate tubes (Beckman Coulter, catalog number: 343775)
4. Coverslips (Sigma-Aldrich, catalog number: Z692263)
5. Microscope slides (Fisher Scientific, catalog number: S17466A)
6. Coverslips (Fisher Scientific, catalog number: S175211A)
7. Guanosine-5'-[(α,β)-methylene]triphosphate (GMPCPP) (10 mM) (Jena Bioscience, catalog number: NU-405S)
8. Tubulin (CYTOSKELETON, catalog number: TL238A)
9. Rhodamine tubulin (CYTOSKELETON, catalog number: TL590M)
10. Dithiothreitol (DTT) (Sigma-Aldrich, catalog number: D9779)
11. N-ethylmaleimide (Sigma-Aldrich, catalog number: E3876)
12. Guanosine 5'-triphosphate (GTP) (Sigma-Aldrich, catalog number: G8877)
13. Hydrogen chloride (HCl) (Sigma-Aldrich, catalog number: 295426)
14. Casein (Sigma-Aldrich, catalog number: C7078)
15. ATP (Sigma-Aldrich, catalog number: A6419)

16. Paclitaxel (Taxol) (Sigma-Aldrich, catalog number: T7402)
17. Piperazine-N,N'-bis(2-ethanesulfonic acid) (PIPES) (Sigma-Aldrich, catalog number: P6757)
18. Magnesium chloride (MgCl₂) (Sigma-Aldrich, catalog number: M4880)
19. EGTA (Sigma-Aldrich, catalog number: E3889)
20. Tris base (Sigma-Aldrich, catalog number: T1503)
21. Potassium chloride (KCl) (Sigma-Aldrich: P5405)
22. Protease inhibitors:
 - Aprotinin (Sigma-Aldrich, catalog number: A3428)
 - Benzamidine (Sigma-Aldrich, catalog number: B6505 or 12072)
 - Note: The product "B6505" has been discontinued.*
 - Pepstatin A (Sigma-Aldrich, catalog number: P5318)
 - Phenylmethylsulfonyl fluoride (PMSF) (Sigma-Aldrich, catalog number: P7626)
 - Leupeptin (Sigma-Aldrich, catalog number: L2884)
 - Soybean trypsin inhibitor (SBTI) (Sigma-Aldrich, catalog number: T6522)
 - Tert-Amyl methyl ether (TAME) (Sigma-Aldrich, catalog number: 283096)
23. Catalase (Sigma-Aldrich, catalog number: C1345)
24. Glucose oxidase (Sigma-Aldrich, catalog number: C6766 or G7141)
- Note: The product "C6766" has been discontinued.*
25. Glucose (Sigma-Aldrich, catalog number: D9434)
26. Phosphate buffered saline (PBS) (Sigma-Aldrich, catalog number: P5368)
27. BRB80 buffer solution (see Recipes)
28. Buffer L (see Recipes)
29. Anti-fade solution (see Recipes)

Equipment

1. Micropipette
 - 2-20 µl (Sigma-Aldrich, catalog number: Z717304)
 - 10-100 µl (Sigma-Aldrich, catalog number: Z717312)
2. Ultra-centrifuge (Beckman Coulter, model: TLA-100)
3. NanoDrop (Thermo Fisher Scientific, Thermo Scientific™, model: NanoDrop™ 2000)
4. Fluorescence microscope (Nikon Instruments, model: Eclipse E600)
5. Water bath (Thermo Fisher Scientific, model: Precision™ General Purpose Baths, catalog number: TSGP02)

Software

1. ImageJ (Version 1.51J, <https://imagej.nih.gov/ij/>)

Procedure

A. Polarity marked microtubules

1. Make bright GMPCPP seed, which serves as a 'primer' for microtubule polymerization as well as an indicator of the microtubule minus-end.
 - a. Make a 10 mg/ml solution of tubulin by dissolving 250 μ g of tubulin into 25 μ l of BRB80.
 - b. Make a 10 mg/ml solution of rhodamine tubulin by dissolving 20 μ g into 2 μ l of BRB80.
 - c. Mix 4 μ l of the 10 mg/ml tubulin solution, 2 μ l of the 10 mg/ml rhodamine tubulin solution, 18 μ l of BRB80, 3 μ l of 10 mM DTT, and 3 μ l 10 mM GMPCPP. Store mix on ice.
 - d. Centrifuge at 410,000 $\times g$ for 7 min at 4 $^{\circ}$ C.
 - e. Take out 3 μ l of bright seed mix supernatant and incubate in water bath at 37 $^{\circ}$ C for 15 min.
2. Make N-ethylmaleimide tubulin
 - a. Make a 10 mg/ml tubulin solution on ice (see step A1a).
 - b. Make a fresh 1 mM N-ethylmaleimide (NEM) solution.
Add 0.5 μ l 50 mM NEM to 25 μ l 10 mg/ml tubulin. Let the solution sit at room temperature for 5 min.
 - c. Add 0.5 μ l of 1 M DTT to 25 μ l NEM-tubulin solution. Let the solution sit on ice for 30 min.
 - d. Centrifuge NEM-tubulin solution at 350,000 $\times g$ for 15 min at 4 $^{\circ}$ C. Collect the supernatant.
 - e. Measure the concentration of supernatant with NanoDrop at 280 nm absorbance.
3. Make dim elongation mix
 - a. To make 10 μ l elongation mix:
 - 1.34 μ l unlabeled tubulin (10 mg/ml)
 - 0.17 μ l rhodamine tubulin (10 mg/ml)
 - 3 μ l NEM-tubulin (4 mg/ml)
 - 1 μ l 10 mM DTT (diluted 1 M DTT into 10 mM DTT first)
 - 0.2 μ l 50 mM GTP
 - 4.29 μ l 1x BRB80
 - b. Incubate the mix in water bath at 37 $^{\circ}$ C for 1 min.
4. Make polarity marked microtubules.
 - a. Add 2 μ l of bright seed mix to 10 μ l dim elongation mix and gently mix.
 - b. Incubate the solution in water bath at 37 $^{\circ}$ C for 30 min, and add taxol stepwise to 10 μ M final solution. Store the microtubules at room temperature and protect from light.

B. Gliding assay

1. Make glass chamber
 - a. Acid wash coverslips (22 x 22 mm) with HCl for 24 h and rinse with ddH₂O for 10 min each, 3 times. Air dry the coverslips for 24 h.
 - b. Attach the coverslip to a glass slide (75 x 25 mm) using two pieces of double-sided tape on the top and bottom. There should be open slits on both sides (Figure 1).

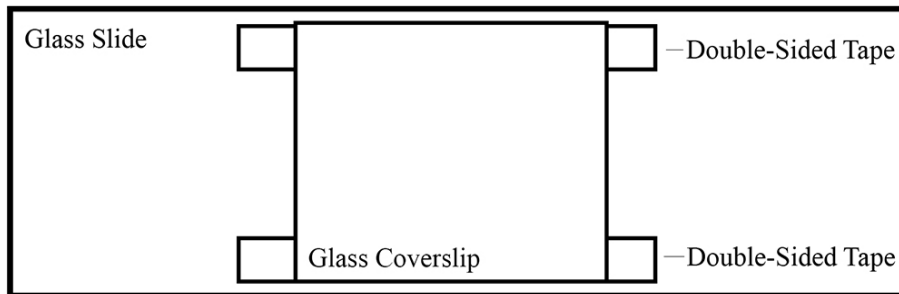


Figure 1. Illustration of glass slide chamber using double-sided tape

- c. Flow in 25 μ l of 1 mg/ml casein into the chamber and let it sit for 5 min. Casein serves as a cushion to anchor motor proteins.
 - d. Wash out the unbinding casein with 3 x 20 μ l buffer L.
 - e. Flow in 10-100 nM purified motor protein (in buffer L), and let it sit for 5 min. Full-length kinesin motor proteins are expressed and purified from baculoviral expression system.
 - f. Wash out unbinding proteins with 3 x 20 μ l buffer L.
 - g. Flow in 25 μ l polarity marked microtubules at a 1/250 dilution into the chamber.
For 100 μ l volume, mix 95.4 μ l buffer L, 1 μ l 100 mM ATP, 1 μ l 1 mM taxol, 2 μ l anti-fade, 0.2 μ l 1 M DTT, and 0.4 μ l MT's.
2. Observe microtubule motility under fluorescence microscope (Nikon E600, excitation wavelength, 547 nm; emission wavelength, 576 nm. Also, see following Video 1).

Video 1. Example of microtubule gliding assay



Data analysis

1. Figure 2 illustrates the principle of a gliding assay. The motor proteins are fixed onto the coverslip and then an aliquot of fluorescent microtubules is introduced into the chamber and caught by the motor. Individual fluorescent microtubules move continuously in a unidirectional

manner. If the motor is plus-end directed, the microtubule itself will translocate in the minus-end leading direction, and vice versa for a minus-end directed motor.

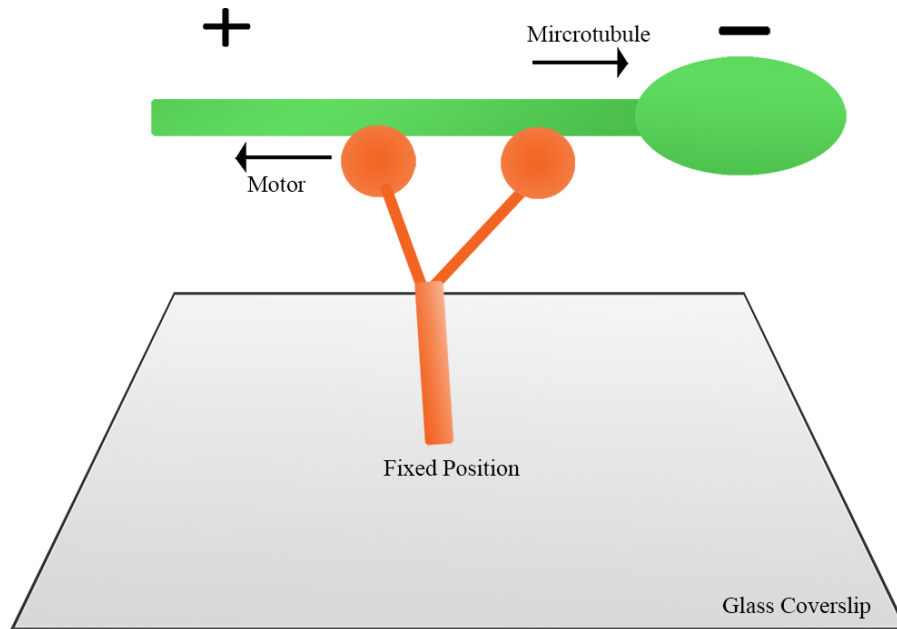


Figure 2. Illustration of a gliding assay. Full-length motor protein is fixed onto coverslip, then polarity-marked microtubules are introduced to the chamber. Microtubule motility is observed under fluorescence microscope.

- Figure 3 shows a plus-end directed motor moving microtubules. The microtubule is in the minus-end leading direction. ImageJ was used to measure the distance of the microtubule movement, and the velocity was subsequently calculated.

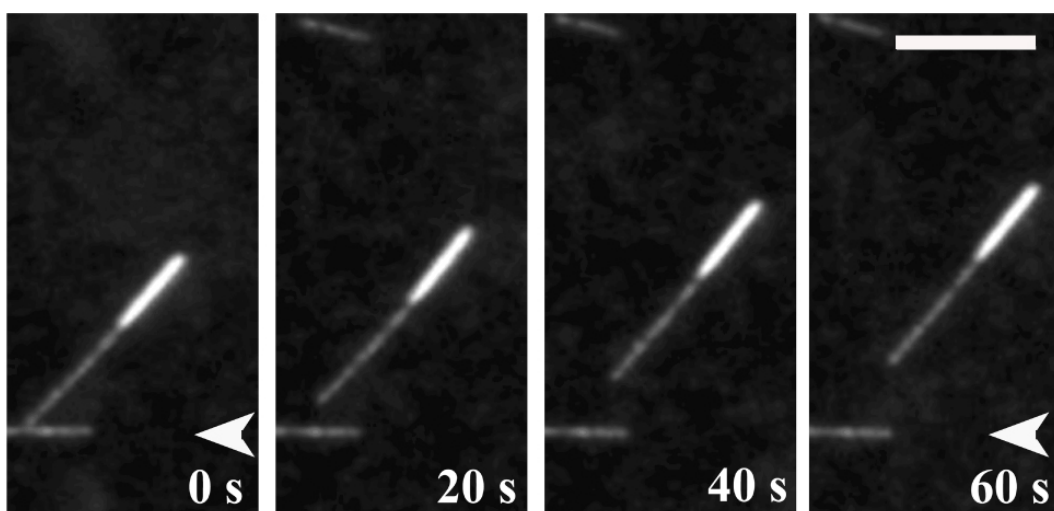


Figure 3. Polarity-marked microtubule translocation in a gliding assay. MT gliding driven by purified kinesin-5 protein. MT is minus-end leading, which suggests that kinesin-5 is a plus-end-directed motor. Scale bar = 5 μm (Tao *et al.*, 2006).

Notes

1. In the gliding assay, the motor proteins must be freshly purified to ensure maximum activity.
2. While making polarity marked microtubules, it is recommended to prepare the bright GMPCPP seed and the dim elongation mix simultaneously. Bright seed mix needs to be added to dim elongation mix right after the seed is made.

Recipes

1. BRB80 buffer solution
 - 1.80 mM PIPES
 - 1 mM MgCl₂
 - 1 mM EGTA, pH 6.8
2. Buffer L
 - 1.20 mM Tris
 - 75 mM KCl
 - 2 Mm MgCl₂
 - 1 mM DTT
 - 0.1 mM ATP
 - Adjust pH to 8.0
 - Protease inhibitors (Final concentrations: Aprotinin, 2 µg/ml; Benzamidine, 20 µg/ml; Leupeptin, 1 µg/ml; Pepstatin A, 1 µg/ml; PMSF, 0.1 mM; SBTI, 100 µg/ml; TAME, 40 µg/ml), added just before use
3. Anti-fade solution
 - a. Mix 10 mg of catalase, 3 mg of glucose oxidase, 180.2 mg glucose, 100 µl 1 M DTT, 1.902 mg EGTA, and 100 µl 10x PBS
 - b. Fill with ddH₂O up to 1 ml
 - c. Aliquot into 20 µl, and store at -80 °C

Acknowledgments

This method is adopted from Tao *et al.* (2016) and Tao *et al.* (2010). The above work was funded by NIH grant GM046409 to W.S. and NIH grant GM 55507 to J.M.S.

References

1. Coppin, C. M., Pierce, D. W., Hsu, L. and Vale, R. D. (1997). [The load dependence of kinesin's mechanical cycle](#). *Proc Natl Acad Sci U S A* 94(16): 8539-8544.

2. Tao, L., Fasulo, B., Warecki, B. and Sullivan, W. (2016). [Tum/RacGAP functions as a switch activating the Pav/kinesin-6 motor.](#) *Nat Commun* 7: 11182.
3. Tao, L., Mogilner, A., Civelekoglu-Scholey, G., Wollman, R., Evans, J., Stahlberg, H. and Scholey, J. M. (2006). [A homotetrameric kinesin-5, KLP61F, bundles microtubules and antagonizes Ncd in motility assays.](#) *Curr Biol* 16(23): 2293-2302.
4. Tao, L. and Scholey, J. M. (2010). [Purification and assay of mitotic motors.](#) *Methods* 51(2): 233-241.